

INTRAcel

FRM™ Fractional RF Microneedling

RF Multi Layer Therapy Outcomes

STRIAE



Photos courtesy of Greg Chernoff, MD

ICE PICK SCARS



Photos courtesy of Jeisys Medical

SKIN TIGHTENING & LIFTING



Photos courtesy of J.C. Poon, MD



Photos courtesy of Marge Uibu, MD

Technical Specification

Applied Energy	RF (Radio Frequency) - 1Mhz
Total RF Energy	Up to 78.4 J/cm2
Spot Size	1 cm x 1 cm
Disposable Tip Matrix	49 electrode needles
Trix Energy Array	Up to 1600 mJ MTZ

System Specification

Display	Color LCD Touch Screen
Weight	138 lb. / 63Kg.
Dimensions	13.7 in. x 15.9 in. x 42.5 in.
Electrical Requirements	AC 110 V, 50~60 Hz, 1 KV



Exclusive distributor of Jeisys Medical in the U.S.A.

Jeisys

Perigee Medical - 844-PERIGEE - info@perigeemedical.com - www.perigeemedical.com - INTC 019816 - © Perigee Medical, All rights reserved.

INTRAcel

FRM™ Fractional RF Microneedling




Jeisys | DISTRIBUTED BY Perigee
medical & aesthetic solutions

The Secret of my Skin

MULTI LAYER THERAPY = HAPPIER PATIENTS

INTRAcel **Multi Layer Therapy** allows precise and independent control of RF energy delivery, depths and thermal zone in any layer of the skin. Physicians can finely tune treatments by selecting and combining the ideal energy needed to achieve the best result:

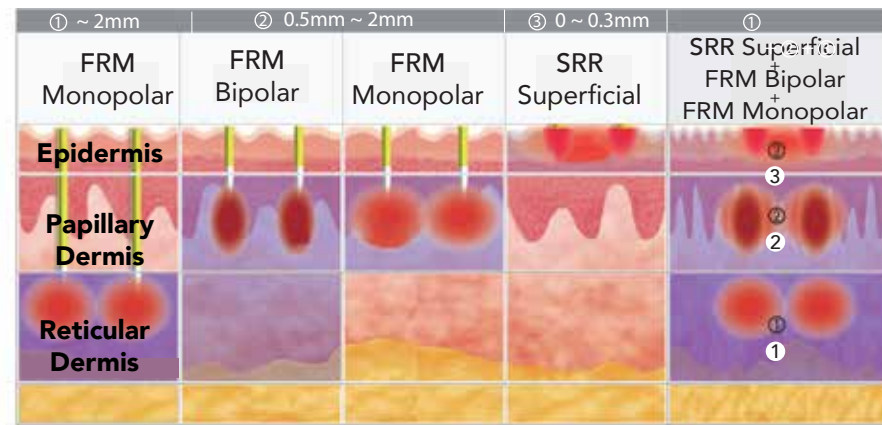


4 modes selection/combination for optimal outcomes

FRM MONOPOLAR for the deeper layers
 +
 FRM BIPOLAR for the papillary dermis
 +
 SRR - MODE (Superficial RF Rejuvenation) for superficial treatments
 +
 T-MODE for conjunction with topicals

INTRAcel BENEFITS

Consistent, more predictable results:

- Precisely deliver energy to the target area
- Adjustable coagulation zone for accuracy and efficacy
- Controllable depth and intensity of treatment
- Active impedance feedback for optimal RF energy delivery
- Greater epidermal protection, reduced risk of PIH
- Minimal downtime for a high degree of patient satisfaction
- Safe on all skin types with minimal pain

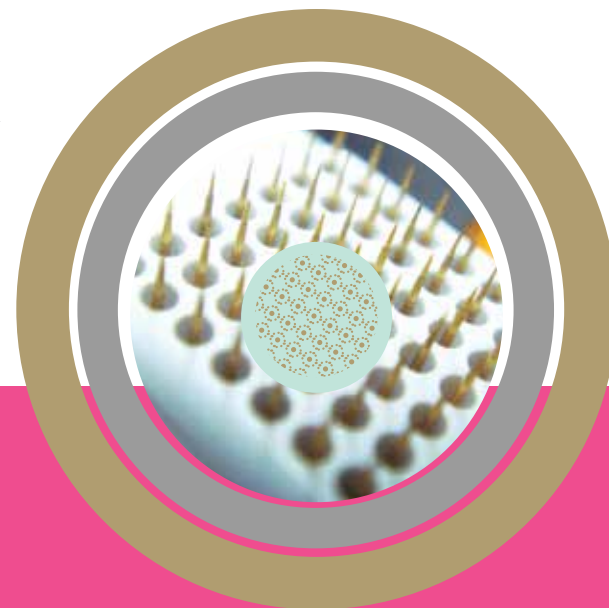


Fractional RF Microneedling™


OPTIMIZED TIPS

(based on 6 years of experience)

INTRAcel tips are optimized to treat various depths of the skin. Each tip has 49 pins specifically designed to reach a selectable depth in the dermis, while the SRR™ sublative tip is optimized to treat the more superficial layers of the skin.

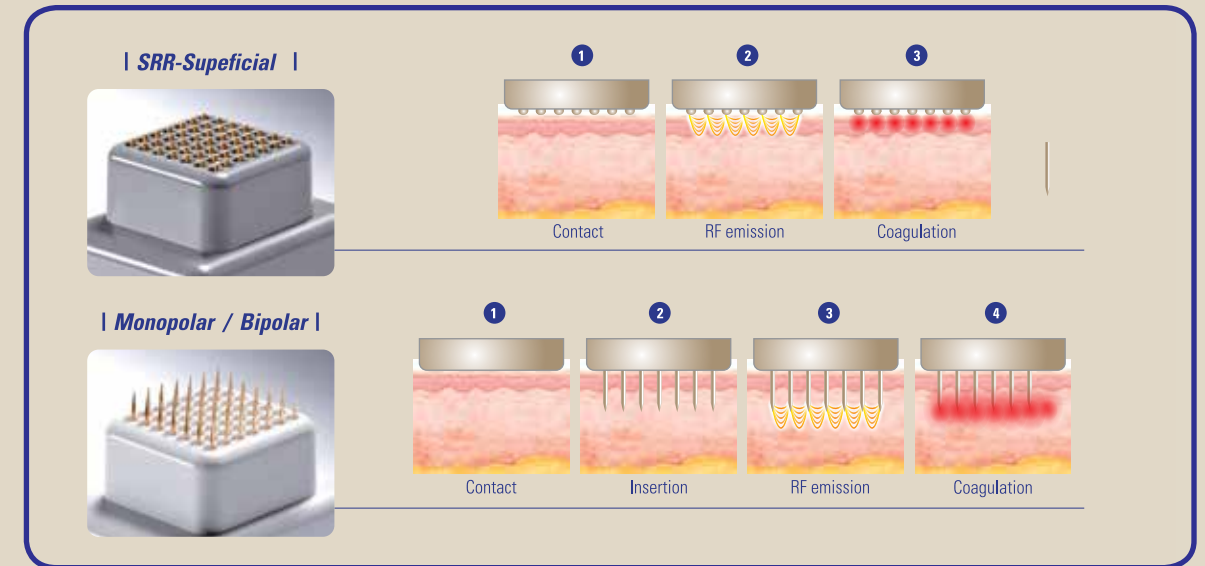

The microneedles are insulated except 0.3 mm at the tip where the RF energy is emitted for maximum patient comfort and targeted energy delivery.

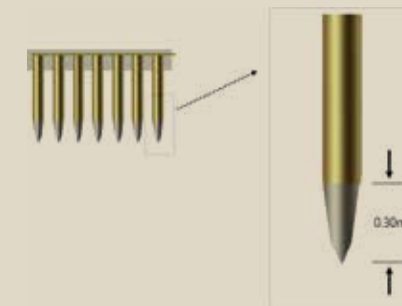

Insulated RF microneedle technology is proprietary technology protected by U.S. Patent 8,666,487.

*"INTRAcel delivers a fractional RF treatment which bypasses the dermal-epidermal junction so we easily minimize risk of PIH and use it on all skin types."**
 - Michael H. Gold, M.D.